

FecalSwab[®]



LIQUID BASED MICROBIOLOGY

FecalSwab[™]

COLLECTION & PRESERVATION
SYSTEM FOR ENTERIC DIAGNOSTICS



Innovating Together[™]

copangroup.com



FecalSwab™ is a **Collection, Preservation** and **Transport** system for Feces and Rectal samples.

It is specifically designed for **enteric pathogens**.



KEY FEATURES

> Ready-to-use

> Sample **preservation**: 48h at RT and 72h at 4°C**

> **Compact container, neat, easy** to transport, use and store

> **Versatile**: sample collected from both defecated stool and rectal swab

> **Validated** for: Culture

> **Supported** analyses: Microscopy *, rapid antigen*, PCR*


> Tested with **molecular platforms***** : Biofire Film Array, BD Max, R Bio-Pharm, GenomEra

> Tube with **maximum sample quantity** indication

> FLOQSwabs® with molded **breaking point**

> Tube with **capture cap**

> **New package**: Vi-pack

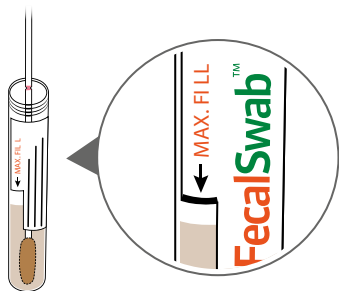

* Analyses applicable after user's validation


** For *C. difficile* 24h at RT and 48h at 4°C

*** Publications available upon request.


SAMPLE WORKFLOW

Microscopy°

Bacteria

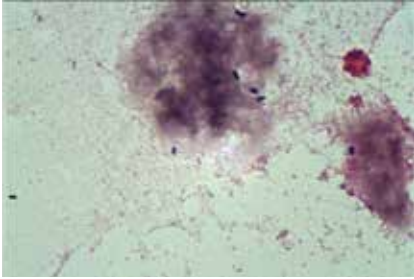

Parasites

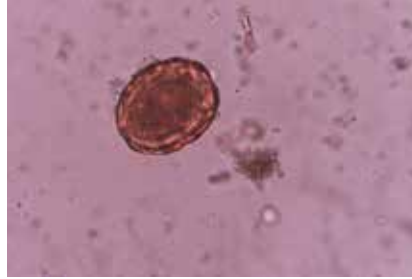

Culture

E.coli
Escherichia coli O157:H
Shigella flexneri
S. typhimurium
E.faecalis (VRE)

V.parahaemoliticus
Y.enterocolitica
C.jejuni
C.difficile



Tested at RT and 4°C

Manual & Automation



Antigen detection°

Adenovirus - *Rotavirus* - *C. difficile*



Molecular assay°

Nucleic acid extraction° & PCR°

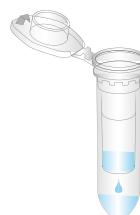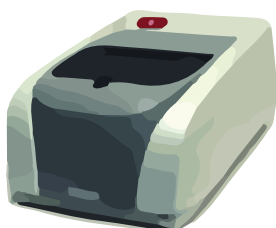



FecalSwab®

ORDERING INFORMATION

Code	Product Description	Pack Size
470CE	FecalSwab™ for manual use, 12x80mm tube filled with 2 ml Modified Cary Blair medium + 1 regular FLOQSwabs®	500 pieces per carton (10 boxes of 50 pieces)
470CE.A	FecalSwab™ for automation, 12x80mm tube filled with 2 ml Modified Cary Blair medium + 1 regular FLOQSwabs®	500 pieces per carton (10 boxes of 50 pieces)
4E020S.A	FecalSwab™ for automation, 12x80mm rounded bottom tube, filled with 2 ml Modified Cary Blair medium + 1 regular FLOQSwabs®	500 pieces per carton (10 boxes of 50 pieces)
4E048S	FecalSwab™ for manual use, 12x80mm tube filled with 2 ml Modified Cary Blair medium + 1 rectal FLOQSwabs®	500 pieces per carton (10 boxes of 50 pieces)
4U031S	FecalSwab™ in bulk for stool specimen: FLOQSwabs® stool transfer devices in resealable pack + 12x80mm tube filled with 2 ml Modified Cary Blair medium in bulk	300 pieces per carton (6 boxes with 50 pieces of tubes + 1 resealable pack of 50 stool transfer devices)



- › European Patent #1608268
- › US Patent #8,114,027
- › US Patent #8,317,728
- › Canadian Patent #2,515,205
- › Japanese Patent #4579902
- › Chinese Patent #101103931
- › Australian Patent #2004226798
- › New Zealand Patent #541560

This document contains information on products which is targeted to a wide range of audiences and may contain product details or information otherwise not accessible or valid in your country. Please be aware that we do not take any responsibility for accessing such information which may not comply with any valid legal process, regulation, registration or usage in the country of your origin.



COPAN ITALIA SPA
via F. Perotti 10
25125 - Brescia, Italy

t. +39 030 2687211
info@copangroup.com
www.copangroup.com

JMKC012B Rev.02 Date 2019/04